

#niwop
f @ in ▶
video.wh.com

Now at
your dealer
or **wh.com**



The new freedom

of rotary
polishing



proxeo ^{TWIST} **Cordless
Polishing System**

Allows you complete **freedom**

Cordless polishing means polishing entirely without restrictions. The system, consisting of disposable prophy angles, handpiece drive and cordless foot control, allows you to change your treatment position in an instant without having to worry about annoying cables. The handpiece is lighter than many smartphones and its ergonomic design adapts perfectly to your movements. This allows you to polish without muscle strain, reach almost every corner of the mouth and always have the best view of the treatment site.



More adaptable than ever

The cordless, ergonomic design means that the handpiece can easily be rotated 360 degrees in the hand, enabling dynamic adjustment of the working position.



Supports your back

Without cable you can pay greater attention to maintaining an ergonomic and back-friendly posture. Cord dragging is now a thing of the past.




Follows you every step of the way

With the cordless foot control you can move much more freely and quickly change your position in relation to the patient.



Allows you complete freedom

Twist and turn the handpiece at will and change your treatment position as needed for better access and an unrestricted view.

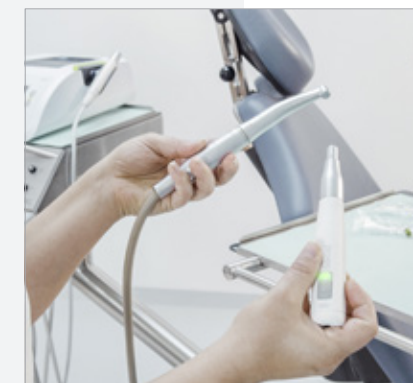


Make it easier on yourself

The cordless handpiece along with the disposable prophy angle is less than half the weight of many connected motors with a speed reduction green contra-angle handpiece – for even greater comfort when polishing.



Never lets you down

The foot control enables you to infinitely vary the handpiece speed and adjust it to the clinical requirements. The durable lithium-ion batteries always ensure sufficient power.

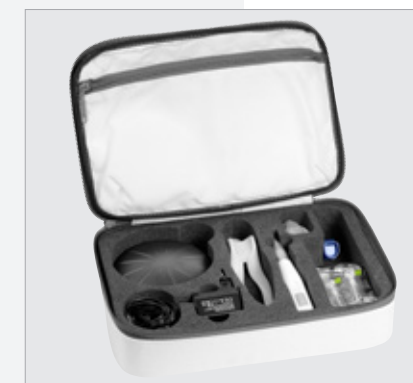

Always there for you

The Proxeo Twist Cordless Polishing System gives you complete freedom when it comes to the treatment location. Simply stow it away in the practical carrying bag and take it with you to the next treatment.

Inspired by Prophy professionals

W&H offers the appropriate disposable prophy angles in various hardness levels for the Proxeo Twist Cordless Polishing System. They were designed in consultation with Prophy users. Great emphasis was placed on features such as optimal adaptation, simple paste pickup and distribution as well as gentle cleaning into the sulcus.




Smooth surface on the edge
for gentle cleaning
into the sulcus.



**Nubs on the
outside**
reduce
splattering of the
paste and polish
the interdental
space at the
same time.



Special internal fins
for easy pickup and targeted
application of paste.

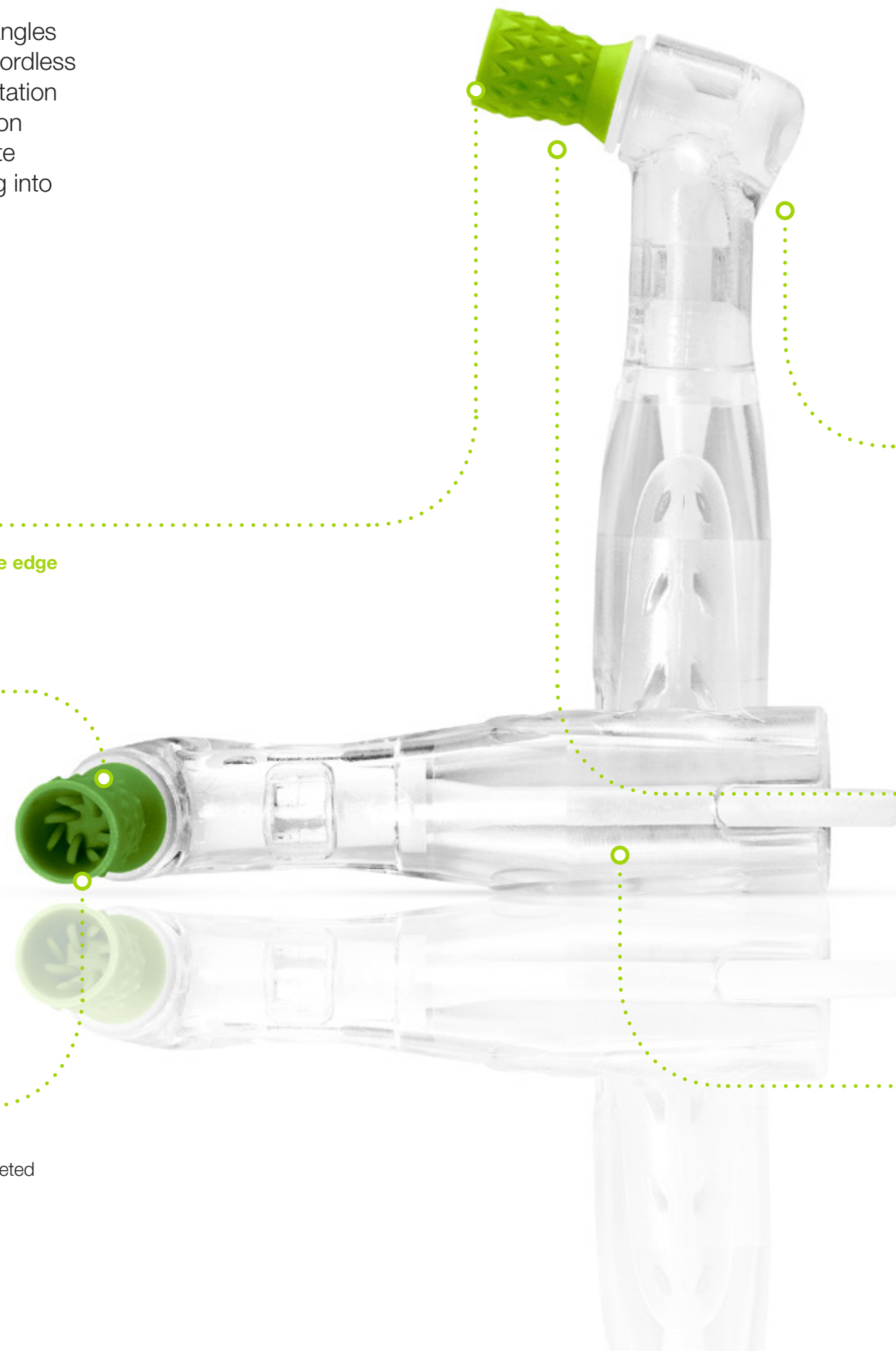

Small head and thin neck
for an optimum view of the
treatment site and better
accessibility of hard-to-reach
areas, such as the molar and
tuber regions.



Optimum adaptation
The cups of the disposable
prophy angles adapt precisely
to the tooth surface.



Makes a good impression
The individually packed clear
cups are opened directly in front
of the patient.



The new lightness during rotary polishing

It's so easy, so easy ...

The new cordless and highly ergonomic Proxeo Twist Cordless handpiece with cordless foot control and the innovative W&H disposable prophy angles for relaxed prophylactic treatments all day long.



Practical holder

You can simply put the handpiece down here between treatments and for recharging.



Easy to reprocess

The smooth surface ensures maximum hygiene. Simply wipe-disinfect the handpiece. The metal sleeve can be removed and sterilised or thermo washer disinfected.

Technical data

Disposable Prophyl Angle	ⓧ Not for reuse	LATEX Latex-free
Angle:	105°	
Connection:	Doriot	
Hardness levels:	soft/apple green firm/dark green	

CE
0297

Cordless handpiece	PL-40 H
Weight:	118 g
Maximum speed:	3,000 rpm
Maximum torque:	2 Ncm
Dimensions (WxDxH):	160 x 25 x 28 mm
Battery type:	Li-Ion
Running time when fully charged:	8 polishing treatments of 6 minutes each
Standby:	automatically after 4 minutes
Charging time:	for full charge (100%) approx. 2 hours for recharge (50%) approx. 30 min

Removable sleeve:

Thermo washer disinfectable Sterilisable

wireless foot control	C-NW
Weight:	190 g
Dimensions (WxDxH):	117 x 117 x 38 mm
Battery type:	Li-Ion
Running time when fully charged:	2 months
Standby:	automatically if not actuated
Charging time:	approx. 3 hours



Bernhard Eder

Product manager W&H, Bürmoos


»The cordless Proxeo Twist offers more freedom of movement and comfort during treatment. We also paid particular attention to the weight. Compared to a contra-angle handpiece with electric motor, the cordless Proxeo Twist weighs less than half as much and therefore minimises wrist fatigue when polishing.«



Monika Riedl

Dental hygienist, Munich

»The cordless Proxeo Twist eliminates the unpleasant pain due to stress and strain usually caused by the tensile strain exerted by a cable.«

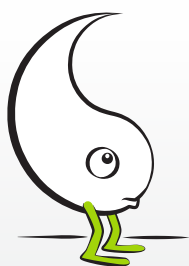

Implantology
without Periodontology is
like Yin without Yang.



#NIWOP: The systematic and evidence-based workflow

W&H Dentalwerk Bürmoos GmbH

Ignaz-Glaser-Straße 53, Postfach 1

5111 Bürmoos, Austria

t +43 6274 6236-0

f +43 6274 6236-55

office@wh.com

wh.com