

#whdentalwerk



video.wh.com

Now at
your dealer
or **wh.com**



Powerful cleaning & disinfection



teon
teon +

Teon

Next step for reprocessing

The new Teon is a highly efficient thermal washer disinfectant. Thanks to its dedicated programs it offers the best cleaning and disinfection for your dental practice. Teon improves the job organization of your hygiene room, reduces the reprocessing time and increases the protection against cross contamination.



Teon is available without any interior, so you can tailor the device exactly to your needs. For example, with the »**starter pack**« (for more information, see page 6).

teon

Teon +

The best approach for reprocessing

Teon +, the premium W&H thermal washer disinfectant, ensures the most demanding clinics an excellent automatic cleaning system for instruments, including handpieces. Teon + offers full traceability to meet current regulatory requirements. Its innovative active air-drying system offers a new standard in terms of quality and performance.



Complete your Teon + with the interior you really need. The device is sold empty, but can be equipped with, e.g. the »**injection pack**« (for more information, see page 6).

teon +

High standard for cleaning and reprocessing



High chamber capacity

Ease of use display

Large selection of accessories

Cleaning of the hollow body

High cleaning performance

Gentle on material thanks to dedicated programs

Traceability and connectivity through Ethernet module

Reliability thanks to heater elements outside the wash and gapless chamber

teon +

Teon + offers in addition:

- › Cleaning and drying of handpieces
- › Active air-drying with HEPA filter




Smooth surface and modern design

The thermal washer disinfectors are distinguished by their high-quality performance and thus ensure a successful, holistic cleaning process for dental instruments. Thanks to their modern, elegant design, they can be perfectly integrated into your reprocessing room.



Safety

The integrated spray arm monitoring enables excellent reprocessing results. The speed variable pump increases/decreases the pressure according to the step of the reprocessing process. In addition, the filter system removes dirt particles to prevent damage occurring to your instruments.

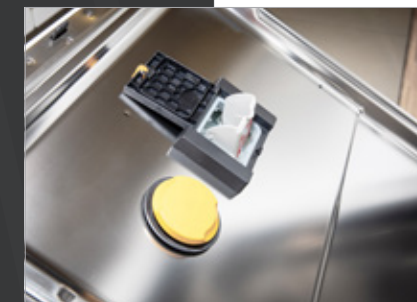

Convenience and ergonomics

The salt container integrated into the door ensures easy and convenient refilling. Thanks to the ergonomics concept, no equipment needs to be removed.



Hygienic rinsing room

The rinsing chamber enables higher throughput with excellent hygiene purity. In addition, the optimal size allows a higher capacity per batch.



Ease of use

The modern design and intuitive user guidance make it particularly easy to operate in the demanding everyday reprocessing routine. The control panel is easy to operate and also serves as a door opener.

Complete and versatile

Reprocessing can be as simple and efficient: with the optimum equipment accessories.

1 – Upper basket

Flexible for inserts and height adjustable.



2 – Insert

Intended for different instruments and for flexible use.



3 – Lower basket

For holding mesh trays as well as various inserts.



4 – Insert

For accommodating various instruments and insertable for upper and lower basket.



»starter pack« 1–6
includes the perfect range
of basic equipment.

»injection pack« 3–8
complemented by an
additional equipment perfectly
and particularly suitable for
handpieces.

5 – Insert

Different compartments
for optimum
loading
of upright
instruments.



6 – Insert

Equipped
with different
compartments
for optimum
loading.



7–8
Upper basket –
Injector: nozzles
For various
instruments with
injector nozzles.



Please note: The devices are
sold empty. Add the appropriate
equipment to your device.



Technical data

CE
0297

Type:	Teon	Teon +
Power supply:	400 V AC; 50–60 Hz/230 V AC; 50–60 Hz	
Current:	16 A	
Overall dimension (w × h × d):	600 × 835 × 600 mm	
Space required for installation (w × h × d):	598 × 835 × 1.200 mm (in total with opened-door)	
Weight (empty):	74 kg	78 kg
Smart Dry:	•	–
Dry Plus:	–	•
HEPA filter:	–	•
Cold water connection:	•	•
Hot water connection:	–	•
Demineralized water connection:	•	•
Auto-close – automatic door lock:	•	•
No. of integrated dispenser pumps for liquid products:	2	1
In-door dosing device for rinse:	•	•
Salt tank – front door access:	•	•
Number of program:	5	6
No. of customizable program:	2	2
Cycle time*:	49'	74' including active air-drying
* time for Vario TD Dental program		

The thermal washer disinfectors Teon and Teon + were designed, certified and validated with the most stringent directives and standards:

93/42/EEC Medical Device Directive	EN 61010-2-040 Safety requirements	DIN EN ISO 15883-1/2 Washer-Disinfector
EN 61326-1 Electromagnetic compatibility	–	EN 1717 Protection against pollution of potable water

W&H Sterilization Srl.
Via Bolgara, 2
24060 Brusaporto (BG), Italy
t +39 035-66 63 000
f +39 035-50 96 988
office.sterilization@wh.com
wh.com

Pictures are for illustrative purposes only. Additional equipment and accessories shown are not included as standard.

20492 AEN Rev. 002 / 21.05.2019
Subject to alterations