



novodent

kuraray
Noritake

PRIRODNA LJEPOTA ZA VAŠE PROTETSKE RADOVE

Saznajte više o KATANA™ ZIRCONIA porodici



KATANA™ Zirconia HTML Plus

High Translucent MultiLayered Plus

Nasljednik HTML-a visoke čvrstoće uz još veću prozirnost

KATANA Zirconia HTML PLUS karakterizira uporaba jedinstvenih materijala koje odgovaraju velikoj brzini sinterizacije i poboljšanja prozirnosti bez ugrožavanja njegove čvrstoće. Napravljen je od kompletno novih materijala i sirovina i koristi Multylayered tehnologiju (četveroslojna struktura). Izvanredan je materijal pogodan za okvire, moderne dizajne, ali i za komplikirane strukture implantata gdje je dosljedna snaga od iznimne važnosti.



KATANA™ Zirconia STML

Super Translucent MultiLayered

Idealan odnos između translucenost i jačine.

Da li mogu cirkonski protetski radovi u punoj formi izgledati prirodno? Da!

S KATANA Zirconia STML blokovima svjetlo se prenosi u incizalno područje i zaustavlja u cervicalnom području. Na ovaj način imate fleksibilnost izrade lijepih protetskih radova u punoj formi, bez obzira na boju nadogradnje. Ovaj proizvod je prvi na tržištu sa karakteristikom cirkona s bojom i prirodnom translucencijom i gradnjom boje. Rezultat rada s ovim diskovima su izrazito estetski i čvrsti protetski radovi.



KATANA™ Zirconia UTML

Ultra Translucent MultiLayered

Translucentnost slična caklini

Radite li ljuškice prirodnog izgleda i ostale protetske radove u prednjoj zoni od cirkona? Da!

KATANA Zirconia UTML ima izrazitu transparentnost, zahvaljujući čemu se postiže prirodna boja zuba te rezultat prirodnostkoju želimo postići. Ova cirkonima karakteristike najveće transparentnosti na tržištu uz paletu raspoloživih prirodnih boja zuba. Proizvod je izvrstan i jedinstven!



KATANA™ Zirconia YML

Yttria MultiLayered

Visoka estetika uz zavidnu čvrstoću

Najnoviji višeslojni disk s dobro uravnoteženom kombinacijom boje/prozirnosti i gradnje snage pri savijanju. Zahvaljujući multi layer tehnologiji razvijen je uskladivanjem tri sirovine različitih udjela itrija. Svi slojevi međusobno savršeno uskladišeni, te su izvršeni svi temeljni koraci proizvodnje, kao što su pročišćivanje, oplemenjivanje cirkonijeve sirovine i dodavanje esencijalnih komponenata poput itrija.

BOJE

HTML Plus	NW	A1	A2	A3	A3.5	A4	B1	B2	B3
	C1	C2	C3	D2	D3				


STML	NW	A1	A2	A3	A3.5	A4	B1	B2	B3
	C1	C2	C3	D2	D3				

UTML	A1	A2	A3	A3.5	A4	B1	B2	B3	B4
	C1	C2	C3	C4	D2	D3	D4		

ENV	EA1	EA2	EA3						
------------	-----	-----	-----	--	--	--	--	--	--

YML	A1	A2	A3	A3.5	A4	B1	B2	B3	
	C1	C2	C3	D2	D3	NV			

30-75% jači od e.maxa
Test savijanja u 3 točke

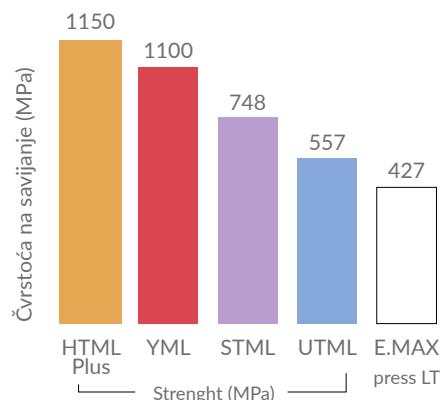

KATANA Zirconia
Multi-Layered serija ima
tri važne karakteristike:
ekonomičnost, brzina i
preciznost. Digitalna
tehnologija daje rješenje
i u estetskoj zoni.

←	14 mm	18 mm	22 mm
	1.380 Kn + PDV 5%	1.740 Kn + PDV 5%	1.960 Kn + PDV 5%

CZR™ FC PASTE STAIN

Inovativni materijal za jednostavnu karakterizaciju restauracija od punog cirkona.

Dostupno u 27 boja, velika mogućnost reprodukcije boja. Široka paleta plavih i sivih boja daje izvrsnu mogućnost transparentnog efekta u incizalnom području. Dva oblika glazure (Glaze i Clear Glaze) olakšavaju mogućnost transparentnosti kroz cijelu krunicu.



Enamel	Dentin	Cervical	Incisal Effect	Body Effect	Pink Tissue
<input type="checkbox"/> Glaze	<input checked="" type="checkbox"/> A+	<input checked="" type="checkbox"/> Cervical 1	<input checked="" type="checkbox"/> Light Gray	<input checked="" type="checkbox"/> Mamelon Orange 1	<input checked="" type="checkbox"/> Salmon Pink
<input type="checkbox"/> Clear Glaze	<input checked="" type="checkbox"/> B+	<input checked="" type="checkbox"/> Cervical 2	<input checked="" type="checkbox"/> Dark Gray	<input checked="" type="checkbox"/> Mamelon Orange 2	<input checked="" type="checkbox"/> Pink
<input type="checkbox"/> Value	<input checked="" type="checkbox"/> C+	<input checked="" type="checkbox"/> Cervical 3	<input checked="" type="checkbox"/> Grayish Blue	<input checked="" type="checkbox"/> Earth Brown	<input checked="" type="checkbox"/> Red
<input type="checkbox"/> Fluor	<input checked="" type="checkbox"/> D+		<input checked="" type="checkbox"/> Blue	<input checked="" type="checkbox"/> Reddish Brown	
			<input checked="" type="checkbox"/> Black	<input checked="" type="checkbox"/> Yellow	
			<input checked="" type="checkbox"/> White	<input checked="" type="checkbox"/> Green 1	
				<input checked="" type="checkbox"/> Green 2	



CZR™ FC Paste Stain Full Kit

Sadrži: Glaze, Clear Glaze (5g cad.); Light Gray, Dark Gray, Grayish Blue, Blue, Cervical1, Cervical2, Cervical3, A+, B+, C+, D+, White, Mamelon Orange1, Mamelon Orange2, Earth Brown, Reddish Brown, Value, Black, Green1, Green2, Yellow, Salmon Pink, Pink, Red, Fluor (3g cad.); ključ boja; ES Liquid (10ml)



CZR™ FC Paste Stain Basic Kit

Sadrži: Clear Glaze (3g.); Dark Grey, Grayish Blue, Cervical2, Earth Brown, A+, White, Salmon Pink (3g .); ključ boja; ES Liquid (10ml)

CZR™ FC Paste Stain Full Kit

(cod. 104-5590-EU)

7.950 Kn
+ PDV 5%

CZR™ FC Paste Stain Basic Kit

(cod. 104-5870-EU)

2.490 Kn
+ PDV 5%

CZR™ Keramika za cirkon



CZR™ Luster CCV Kit

Sadrži: Lt0 (10G), Lt1 (10G), T Blue (10G), Aqua Blue 1 (10G), Aqua Blue 2 (10G), Lt Natural (10G), Lt Super Gray (10G), Creamy Enamel (10G), Sun Bright (10G), Incisal Aureola (10G), Creamy White (10G), Lt Yellow (10G), Ccv-1 (10G), Ccv-2 (10G), Ccv-3 (10G), Ccv-4 (10G). Ccv-2 (10G), Ccv-3 (10G), Ccv-4 (10G).



CZR™ Internal Stain Complete Kit

Sadrži: White (3G), Incisal Blue 1 (3G), Incisal Blue 2 (3G), Mamelon Orange 1 (3G), Mamelon Orange 2 (3G), Cv-1 (3G), Cv-2 (3G), Cv-3 (3G), Earth Brown (3G), Reddish Brown (3G), Fluoro (3G), Salmon Pink (3G), Red (3G), A+ (3G), B+ (3G), C+ (3G), D+ (3G), Gray (3G), Bright (Dilution) (3G), IS Liquid (10 ML), Ključ Boja.



CZR™ Trial Kit SS A2

Sadrži: OBA2 (5g), A2B (5g), E2 (5g), LT1 (5g), IS Liquid (10ml), Forming Liquid (10ml), CZR PRESS SBS SSA2 (3g).

CZR™ Luster CCV Kit

(cod. 104-1790-EU)

2.322 Kn
+ PDV 5%

CZR™ Internal Stain Complete Kit

(cod. 104-1840-EU)

3.600 Kn
+ PDV 5%

CZR™ Trial Kit SS A2

(cod. 104-5530-EU)

300 Kn
+ PDV 5%

EX-3™ Keramika za metal



EX-3™ N Color Intro Kit

Sadrži: Paste Opaque NA2 (6g), Paste Opaque NA3 (6g), Body NA2B (10g), Body NA3B (10g), Enamel E2 (10g), Enamel E3 (10g), Luster LT1 (10g), Luster Creamy Enamel (10g), Luster CCV-1 (10g), Translucent ADT (10g), Translucent TX (10g), Forming Liquid (10ml)

EX-3™ N Color Intro Kit

(cod. 102-0050-EU)

1.429 Kn
+ PDV 5%

EX-3™ NP Bonder

6 g
(cod. 102-4532-EU)

241 Kn
+ PDV 5%

EX-3™

10 g

94 Kn
+ PDV 5%

*10 x snažniji na kobalt-kromu
Adheziv koji poboljšava vezu
između keramike i plemenitih
legura, na bazi krom-kobala

* opcije boja
Opacious Body OBA3
Clear Cervical CCV-2
Speed Enamel S1
Speed Enamel S2
Speed Enamel S3
Speed Enamel S4

PANAVIA Cementi



PANAVIA V5

Panaviu V5 odlikuje iznimna čvrstoća i naglašena estetska komponenta!

Panavia V5 je najjači cement ikad razvijen, a dostupan je u 5 nijansi. Koristi se za cementiranje krunica, mostova, ljuštica izrazito je pouzdan za nagoradnje na implantatima. Postupak započinjete s tretiranjem zuba s PANAVIATM V5 primerom, nastavljate s tretiranjem protetskog rada s CLEARFILTM CERAMIC PRIMER PLUS, potom cementirate - i to je to. PANAVIATM V5 AUTOMIX CEMENT olakšava rukovanje i apliciranje uz jednostavno uklanjanje viška cimenta. Automix šprica omogućava brz, jednostavan i precizan rad - a konzistentno optimiziran cement omogućuje da postavite protetski rad na način koji želite.

PANAVIA™ V5 Professional Kit

Tooth Primer (2 ml), Ceramic Primer Plus (2 ml), Panavia V5 Pasta u 5 boja (2.4 ml/4.2g), probne paste u 5 boja, ukupno 12 ml u 5 boja.

2.790 Kn
+ PDV 5%

PANAVIA™ V5 pasta

pasta 4,6 ml + 20 mješalica

890 Kn
+ PDV 5%

PANAVIA™ V5 Uvodni set

Boja Universal (A2) ili Clear Tooth Primer (2 ml), Ceramic Primer Plus (2 ml), Panavia V5 Pasta (2.4 ml/4.2g)

950 Kn
+ PDV 5%



PANAVIA SA CEMENT UNIVERSAL

Samoadhezivni kompozitni cement

Jedinstveni i inovativni kompozitni cement, samoadhezivni, indiciran za sve vrste materijala (keramiku, litij disilikat, kompozit) bez potrebe silaniziranja. Uz jedinstveni MDP monomer u formuli je dodana komponenta LCS (LONG CHAIN SILANE) koja osigurava snažnu kemijsku vezu na sve vrste površina. To je samojetkajući cement koji ne treba korištenje primera i može se koristiti u svim slučajevima – suhom i vlažnom okruženju; kod cementiranja zirkona, kompozita, litij disilikata (E-max krunica). Zahvaljujući MDP monomeru cement ima funkciju samoadhezije na bazi kemijske veze. U cement je uključena jedna nova komponenta LCSi MONOMER (LONG CHAIN SILANE MONOMER). To je jedinstvena formula koju nema niti jedan drugi cement – što znači da nije potreban silan za krunicu (kao Ceramic Primer ili slično).

PANAVIA™ SA Cement Universal

Boja A1 ili translucent,
1 x tubica 8,2 g, sa priborom

780 Kn
+ PDV 5%

**VALUE
PACK
2+1**

PANAVIA™ SA Cement (3kom)

Boja A1 ili translucent,
1x tubica 8,2 g, sa priborom

1.560 Kn
+ PDV 5%



ALLOY PRIMER

Primer za kondicioniranje svih vrsta metala u protetskim radovima. Jednostavan za korištenje i aplikaciju. Nanosi se na metalnu površinu u tankom sloju i pričeka 5 sekundi prije početka cementiranja.

ALLOY PRIMER 5 ml

510 kn
+ PDV 5%



KATANA CLEANER

Prvi univerzalni intra i ekstra oralni restaurativni cleaner površine.

KATANA™ Cleaner je prvo univerzalno sredstvo za čišćenje svih površina koji ima pH 4.5 što omogućuje ne samo extraoralno već i intraoralno čišćenje.

Indikacije:

- Intraoralna primjena: zub, bataljak, korijenski kanal
- Ekstroronalna primjena: cirkonske strukture, litij disilikat, keramika, materijali od kompozita CAD-CAM krunice, kompozit, metali (legure plemenitih / neplemenitih metala), kolčići (kompozitni i metalni)

KATANA™ Celaner
4 ml
198 Kn
+ PDV 25%



Clearfil Ceramic Primer plus

Novi Clearfil Ceramic Primer se odlično slaže sasvim vrstama keramike. Osim već dokazanog monomera MDP, koji osigurava snažnu vezu, jednokomponentni primer takođe sadrži etanol i silan agent MPS, koji osigurava snažnu vezu na sve materijale od kojih se izrađuju estetske restoracije. Kroz optimalnu kombinaciju ta tri elementa, Clearfil Ceramic Primer održava izvrsnu dugotrajanu adheziju.

Clearfil Ceramic Primer plus
4 ml
690 Kn
+ PDV 25%



INTRAORALNO



EKSTRAORALNO



**DEKONTAMINACIJA
ZA 10 SEKUNDI**

KATANA™ TWIST DIA 3 koraka za blistav sjaj zirkona



KATANA™ TWIST DIA
3 kom
255 Kn
+ PDV 25%



KATANA™ ZIRCONIA BLOCK STML

Multi-layered KATANA STML
Zirconia Block sa 4 razine
cirkona u gradiranim bojama.

Proizvod namijenjen za rad u ordinaciji u prirodnim bojama – efikasniji i brži način rada. Jačina 763 MPa!
15 minuta obrada + 18 minuta sintranje

Dostupno u 3 veličine: 12Z, 14Z i 14Z L
Dostupno u bojama NW, A1, A2, A3, A3.5, A4, B1, B2, B3, C1, C2, C3, D2 i D3.
Dostupno za CEREC.

Indikacije:

- Krunice u prednjoj i stražnjoj zoni
- Ljuskice
- Inlej - onlej
- Most 3 članova



		a (height)	b (width)	c (length)
12Z	Before sintering	15.3 mm	19.2 mm	20.2 mm
	After sintering	12.2 mm	15.4 mm	16.2 mm
14Z	Before sintering	17.8 mm	19.2 mm	20.2 mm
	After sintering	14.2 mm	15.4 mm	16.2 mm
14Z L	Before sintering	17.8 mm	19.2 mm	40.0 mm
	After sintering	14.2 mm	15.4 mm	32.0 mm

12Z (5 kom)
990 Kn
+ PDV 5%

14Z (5 kom)
1.050 Kn
+ PDV 5%

14Z L (5 kom)
1.650 Kn
+ PDV 5%

MULTI
LAYERED

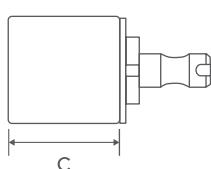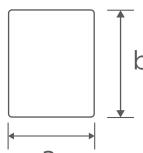

KATANA™ AVENCIATM

Kompozitni blokovi
Hibridna keramika



Dostupno u 2 veličine:
12 (jedna krunica) i 14L (most)
Boje: A2, A3 i A3.5
Dostupno za: CEREC i UNIVERSAL nastavke

Indikacije:

- Inlej - onlej
- Premolarne krunice
- Krunice u prednjoj i stražnjoj zoni
- Most 3 članova

	a (height)	b (width)	c (length)
12	10 mm	12 mm	15 mm
14L	14.5 mm	14.5 mm	18 mm

12 Universal
(5 kom)

790 Kn
+ PDV 5%

14L Universal
(5 kom)



840 Kn
+ PDV 5%

12 Cerec
(5 kom)

870 Kn
+ PDV 5%

14L Cerec
(5 kom)

920 Kn
+ PDV 5%

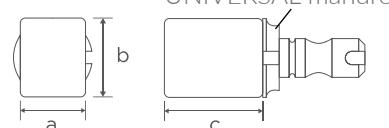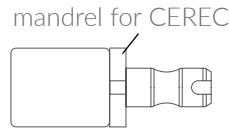

PERFECTA 900

Laboratorijski motor za cirkon
U osnovnoj konfiguraciji je turbina,
a dodatno se može spojiti i motor LA-66
i time dobiti kompletno rješenje za lab

35%

PERFECTA 900
16.400 Kn
10.660 Kn
+ PDV 5%

35%

Motor LA-66
za ugradnjnu na
PERFECTA 900
8.025 Kn
5.216 Kn
+ PDV 25%

35%

PERFECTA 600
15.200 Kn
9.880 Kn
+ PDV 25%

PERFECTA 600

Laboratorijski motor za cirkon
sa motorom LA-66

- motor 7.8 Ncm / 1.000 - 50.000 rpm
- Air-jet funkcija
- 3 verzije: stolna, koljenska i nožna
- 24 mjeseca garancije

PERFECTA 300

Laboratorijski motor

30%

PERFECTA 300
9.800 Kn
već od **6.860 Kn**
+ PDV 25%

Uvoznik & distributer

Novudent d.o.o. • info@novudent.hr • www.novudent.hr

MATULJI Cesta dalmatinskih brigada 7 • 51211 Matulji • Tel: 051 701 645

ZAGREB Zavrtnica 3 • 10000 Zagreb • Tel: 01 457 7846

SPLIT Tel: 091 2 704 712



INFO & NARUDŽBE
091 270 4705